Gebruiksaanleitung Manual Mode d'emploi

(D)

Seite 2

(GB)

Page 4

(F)

Page 6

Alpha[®] Einstieghilfe Ultra
Alpha[®] Transfer aid Ultra
Alpha[®] Aide au Transfert Ultra

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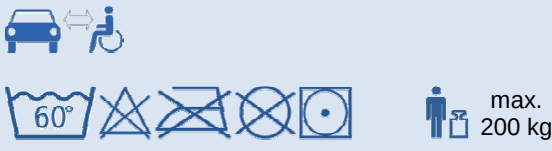
1. Specifications and Scope of Supply

Item no.:	PM-8100
Dimensions:	47 x 45 cm
Material:	Velour/Nylon, latex free
Colour:	grey/red
ILN:	4250355381002
Storage:	-5 – 40° C protect from sunlight
Scope of supply:	
• Transfer aid Ultra • Manual	
Made in Germany	
	
max. 200 kg	
CE	

2. Use / Placement

The outer surfaces of the transfer aid Ultra are made of velour whilst the inner surfaces are made of a woven high-tech material which only permits lateral movement. Forward movement and rearward movement in the car are thus inhibited.

The transfer aid Ultra makes it easier for the user to get in and out of the car without assistance. It also makes it possible for immobilised persons to get in and out painlessly with assistance.

It is placed on the car seat with the rough surface of the smaller seat cushion facing down. The straps are fastened round the backrest. Pull them tight.

This product is suitable for re-use.

3. Maintenance / Care

The product does not require any special maintenance. It can be washed in a washing machine at 60° C and is suitable for low dryer temperatures (45 – 50° C). When using non-standard washing agents, bleaching agents, softener or other chemical washing additives, the functional properties of the product can strongly degrade. The robust material can be disinfected with the usual wipe disinfectants. Check the compatibility before use where required on the basis of the disinfectant's instructions. We recommend the use of alcohol- and aldehyde-free agents, e.g. Trionic®. There is also some care information on the label.

4. Warranty

- Petermann gives a warranty that the product is devoid of material defects and manufacturing defects.
- The warranty is valid for a period of **3 years**. It does not cover any defects which are due to improper use.
- Petermann cannot be held liable for any accidental or indirect injury or damage.

5. Visual and Function Inspection / Safety Information / Re-use

Repositioning products are also subject to natural aging. Check the equipment regularly and carefully to make sure it is working properly and especially before re-use. Make sure that the product is in good condition and that the material and the seams are not damaged. Check the function of the product by applying pressure and loading to the grips and the clasp and make sure that the material resists the load. The smooth surface of the low-friction material makes it easy to slide and move. The high-friction has a rough surface to prevent slippage. The best way to test these properties is by comparing the product with a new product.

For re-using the product it must be washed and disinfected before, see 3.

WARNING:

- ***Faulty products must no longer be used.***
- ***The transfer aid may only be used together with a car safety belt while driving in a vehicle.***
- ***If the patient's mobility is restricted, it is essential to supervise the motions of transfer by providing trusted help.***
- ***The function of the product may be adversely affected if it is washed with non-standard detergents or bleaches, fabric softeners or other chemical detergent additives.***
- ***The smooth surface considerably minimises resistance to friction during the change of position. Carefully check the transfer surfaces, making sure they are stable and brakes are applied (e. g. wheelchair, bed).***

6. Product Labeling

The Petermann Alpha® transfer aid Ultra is a medical product, class I and complies with the MDD directive 93/42/EEC.

Petermann GmbH is certified in accordance with the directive DIN EN ISO 13485 by the German TUV.

7. Disposal Information

For the transfer aid Ultra there are no special disposal rules; it can be disposed of in accordance with national and local regulations as necessary.

If you have any questions about the use of our products, please call our service number:

+49 (0) 9868 9339-0

Petermann GmbH • Ziegelhaus 12 • D-91601 Dombühl
Postal address: Petermann GmbH • Waldweg 3 • D-91601 Dombühl
Phone: +49 (0) 9868 9339-0 • Fax: +49 (0) 9868 9339-20
E-mail: petermann@pm-med.de • www.pm-med.de

"Information subject to technical modifications to our products and deviations in colour or printing errors"

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Petermann GmbH est certifié par le TUV conforme à la directive DIN EN ISO 13485.
Adresse postale: Petermann GmbH • Waldweg 3 • D-91601 Dombühl